


# RENOVO | RX

## Delivering Therapy Where it Matters™




Introducing the RenovoCath® Device Designed  
for Targeted Diagnostic & Therapeutic Agent Delivery

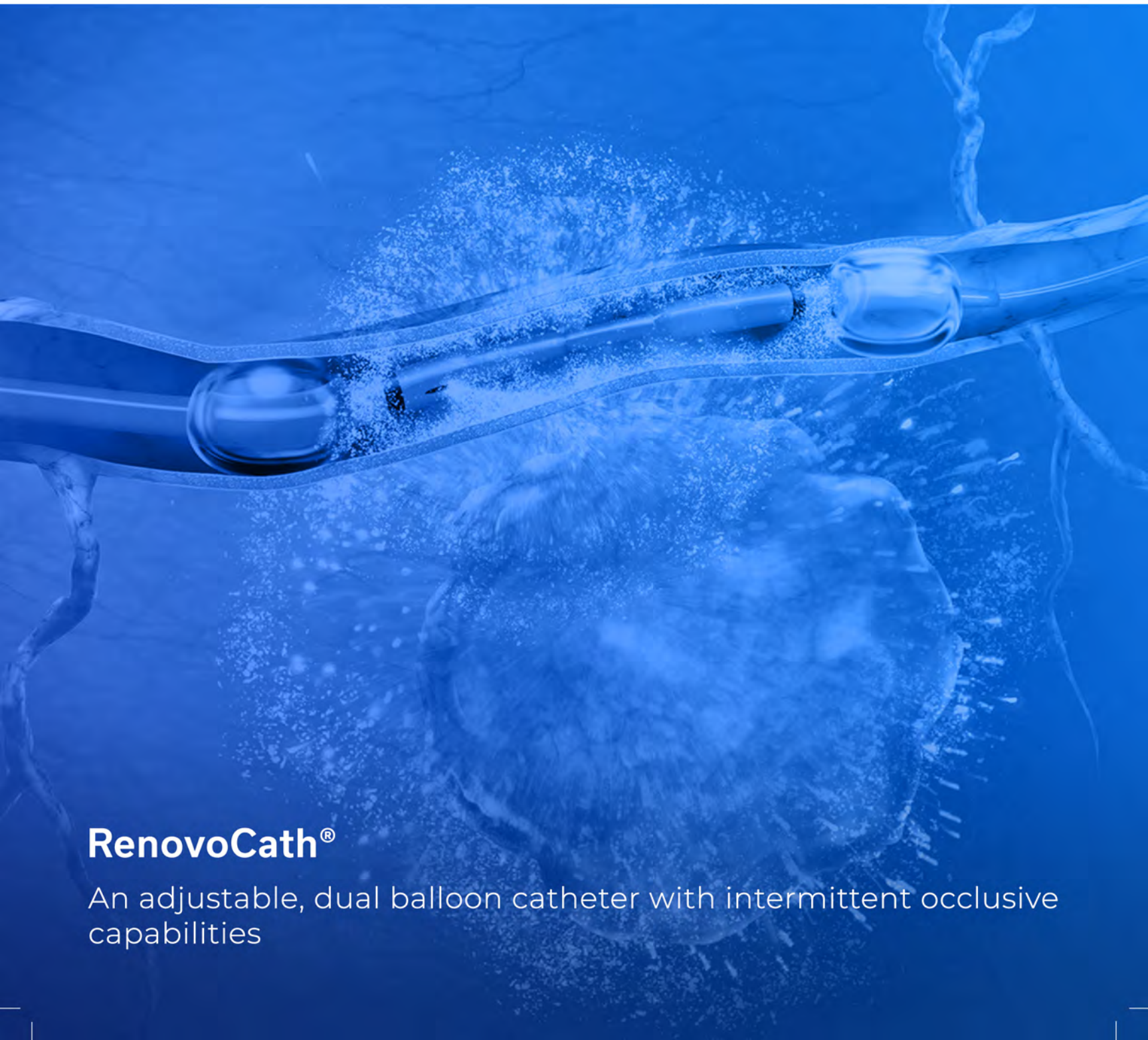

**RenovoCath®**


An adjustable, dual balloon catheter with intermittent occlusive capabilities

**RenovoCath** is a patented, FDA-cleared delivery system enabled by an adjustable, dual balloon catheter with intermittent occlusive capabilities. RenovoCath enables physicians to isolate sections of the blood vessel through the adjustment of the distance between the occlusion balloons, thereby excluding any off-target side branches and/or including tumor feeders.

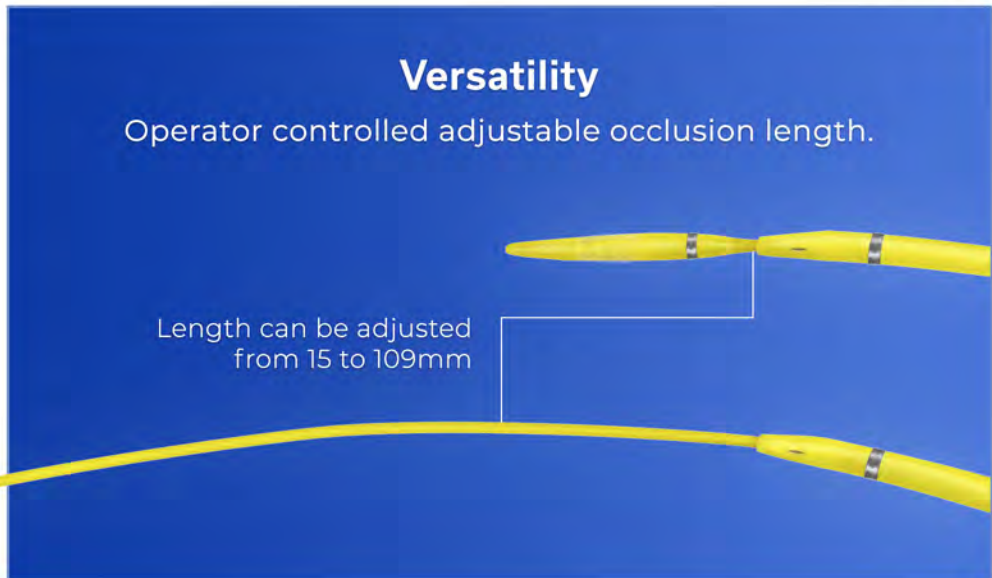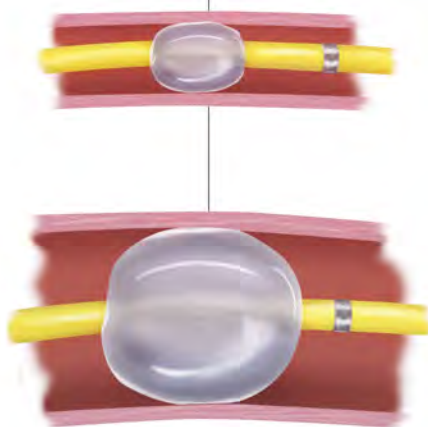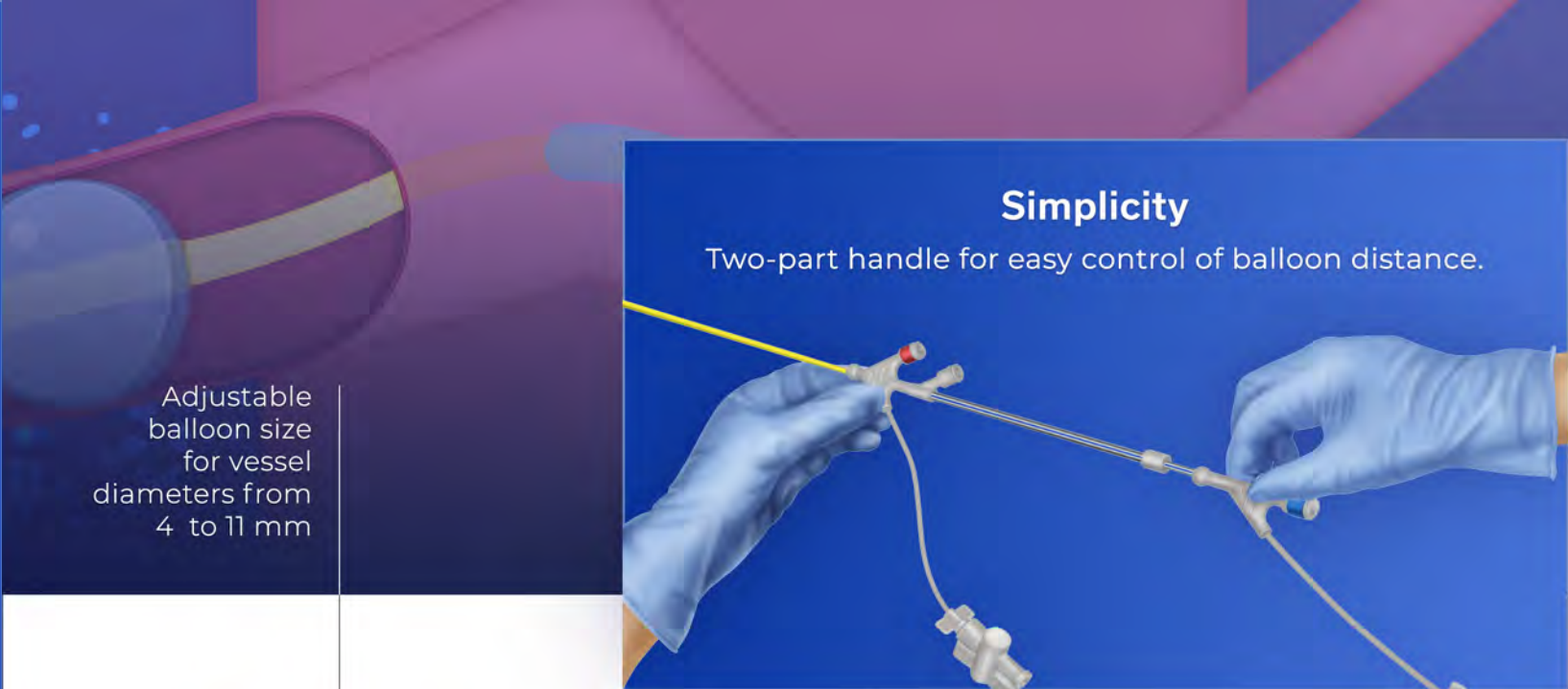

| Balloon Compliance            |                      |
|-------------------------------|----------------------|
| Inflation Volume (cc) Average | Vessel Diameter (mm) |
| 0.10                          | 3                    |
| 0.10                          | 4                    |
| 0.14                          | 5                    |
| 0.19                          | 6                    |
| 0.28                          | 7                    |
| 0.39                          | 8                    |
| 0.56                          | 9                    |
| 0.69                          | 10                   |
| 0.99                          | 11                   |

NOTE: Maximum volume tested is 1.07 cc

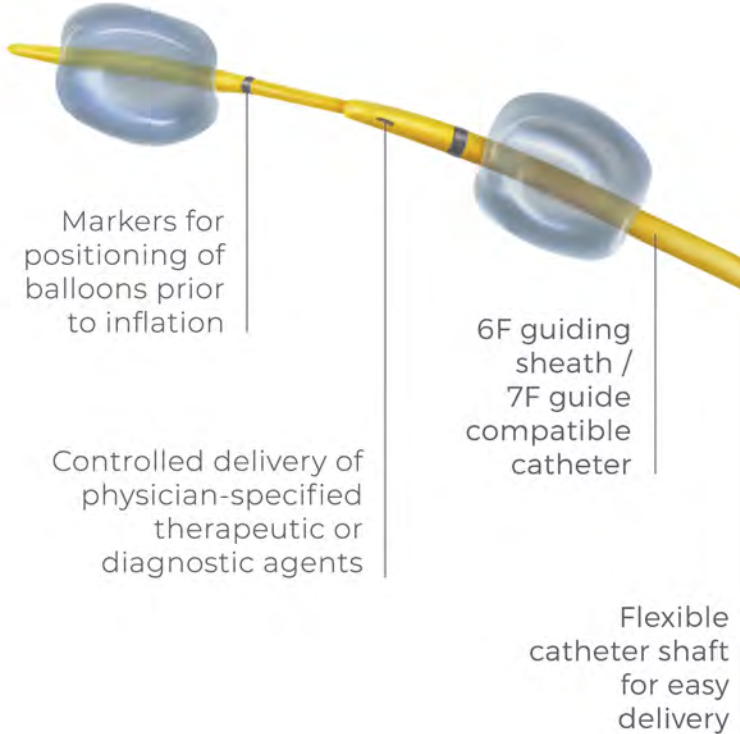

| RenovoCath Delivery System Specifications                                                                                                                                                       |                             |                            |               |                              |                          |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------|----------------------------|---------------|------------------------------|--------------------------|
| Balloon Diameter                                                                                                                                                                                | Balloon to Balloon Distance | Effective Length*          | Catheter O.D. | Guiding Sheath Compatibility | Guide Wire Compatibility |
| 3 – 11 mm (adjustable)                                                                                                                                                                          | 15 – 109 mm (adjustable)    | 75 cm – 85 cm (adjustable) | 0.076 inches  | 6F                           | 0.014 in                 |
| * Effective length is the length of the catheter that can be inserted into the body and is measured from the distal edge of the strain relief to distal tip. The overall length is approximate. |                             |                            |               |                              |                          |

# Delivering Therapy Where it Matters

## Targeted Diagnostic & Therapeutic Agent Delivery

RenovoRx is committed to transforming the lives of patients by providing innovative solutions to enable targeted delivery of diagnostic and therapeutic agents. RenovoCath is FDA-cleared and available for commercial sale.

### RenovoCath Orders & Inquiries

 650.284.4433    650.397.4433    [orders@renovorx.com](mailto:orders@renovorx.com)   Catalogue number – MMRC120-DB-1111

### Indications for Use

RenovoCath is intended for the isolation of blood flow and delivery of fluids, including diagnostic and/or therapeutic agents, to selected sites in the peripheral vascular system. RenovoCath is also indicated for temporary vessel occlusion in applications including arteriography, preoperative occlusion, and chemotherapeutic drug infusion. RenovoCath is intended for general intravascular use in the peripheral vasculature in arteries 3mm and larger. The RenovoCath is intended for use in arteries from 3mm in diameter for vessel entry and to occlude vessels ranging between 3mm to 11mm in diameter.

The diagnostic and/or therapeutic agents are to be used in accordance with specifications outlined by the respective agent manufacturer.

For detailed information on RenovoCath including warnings, precautions, and contraindications, please refer to the product IFU on the RenovoRx website (<https://renovorx.com/for-clinicians/>).

**RENOVO | RX**

 2570 W El Camino Real, Mountain View, California 94040    650.284.4433    [renovorx.com](https://renovorx.com)